

Stainless Dipper

GMP Correct Liquid Sampler

- 316 stainless steel construction
- Fully welded
- Satin finish

| Stock No. | Nominal Volume | Length | Cup Diameter |
|--------------|----------------|--------|--------------|
| 675 5005A100 | 100ml | 600mm | 39mm |
| 675 5005B100 | 100ml | 1000mm | 39mm |
| 675 5005A250 | 250ml | 600mm | 64mm |
| 675 5005B250 | 250ml | 1000mm | 64mm |
| 675 5005A500 | 500ml | 600mm | 90mm |
| 675 5005B500 | 500ml | 1000mm | 90mm |




[QUICK QUOTE LINK](#) [ONLINE PRICING LINK](#)

PTFE Dipper

Excellent Chemical Resistance - GMP correct

- PTFE product contact parts
- Screw-in extension rods to increase the length of the handle

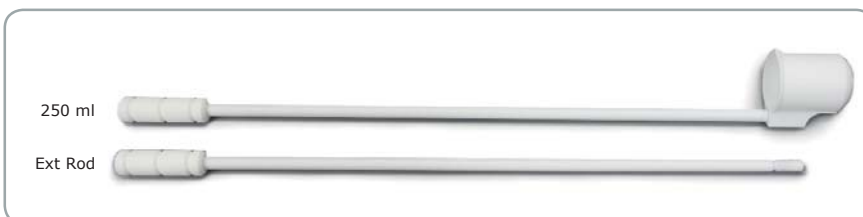

| Stock No. | Nominal Volume | Cup Diameter | Handle |
|---------------|-----------------------------------|--------------|------------------------|
| 675 5100A100 | 100ml | 55mm | 600mm long - removable |
| 675 5100A250 | 250ml | 70mm | 600mm long - removable |
| 675 5100A500 | 500ml | 80mm | 600mm long - removable |
| 675 5100A1000 | 1000ml | 105mm | 600mm long - removable |
| 675 5100E600 | Handle Extension rod - 600mm long | | |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)



Chemical Dipper

Robust Polypropylene Dipper

- Ideal for general industrial applications
- Durable polypropylene construction

| Part No. | Beaker Volume | Overall Length |
|---------------|---------------|----------------|
| 675 5065F250 | 250ml | 1000mm |
| 675 5065F500 | 500ml | 1000mm |
| 675 5065F1000 | 1000ml | 1000mm |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)



Utility Dipper


Simple Liquid Sampler

An all stainless steel sampler suitable for sampling liquids and suspensions.

- Takes a sample at the surface of the liquid
- Made from 304 Stainless steel

Operation

1. Insert the sampler into the product. The liquid will flow into the sampler.
2. Withdraw the sampler from the product
3. Raise the centre control rod to discharge the sample

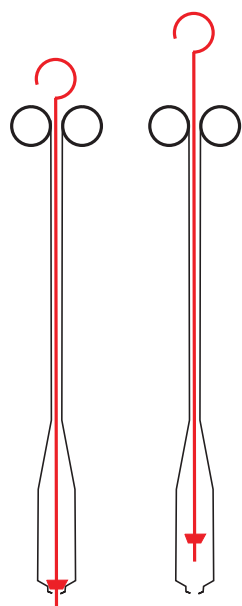

Raise the centre control rod (shown in red) to discharge the product



| Stock No. | Length | Volume | Cup Diameter |
|--------------|--------|--------|--------------|
| 675 5170A50 | 460mm | 50ml | 32mm |
| 675 5170A100 | 540mm | 100ml | 32mm |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)

Multi Dipper 2


Multilevel Liquid Sampler

The Multi Dipper 2 is the ideal solution for taking liquid samples from different levels. The sample collection bottles remain closed until the operator uses the control chain to remove the stoppers.

High security stainless steel shackles are used to secure the plugs to the control chain

Product contact parts are:

- Sampler frame 316ss
- Control Chain 316ss
- Plugs 316ss
- Plug Shackles 316ss
- Sample Bottle HDPE

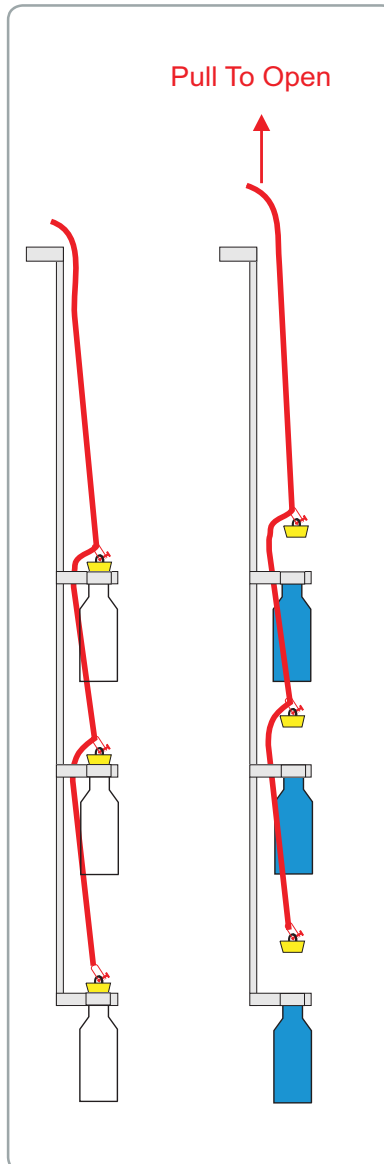

Operation

1. Screw the bottle into the holders
2. Insert the plugs into the tops of the sample bottles
3. Connect the Control Chain to the plugs
4. Insert the assembled sampler into the liquid to be sampled
5. Pull on the Control Chain to open the bottles
6. Remove from the liquid tank

Special sizes and designs are available on request

| Stock No. | Length | Bottles |
|---------------|------------------------------------|-----------------------------------|
| 675 5045A1000 | 1000mm | 3 x 250ml included (675 8165A250) |
| 675 5045A1500 | 1500mm | 3 x 250ml included (675 8165A250) |
| 675 5045A2000 | 2000mm | 3 x 250ml included (675 8165A250) |
| 675 8165A250 | Spare sample bottle - HDPE - 250ml | |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)



Enviro Dipper

The Enviro Dipper is a simple and yet effective method of taking a surface sample. Simply lower the bucket into the liquid.

The stainless steel bucket can be ordered with either a stainless steel chain or a braided polypropylene cord.

| STOCK NO | DESCRIPTION |
|-------------|--|
| 675 A724300 | Stainless Steel Bucket - 300ml - 304ss |
| 675 5070S10 | Chain - 316 stainless steel - 10m - with 316ss shackle |
| 675 5070S25 | Chain - 316 stainless steel - 25m - with 316ss shackle |
| 675 5070S50 | Chain - 316 stainless steel - 50m - with 316ss shackle |
| 675 5070P10 | Braided Cord - Polypropylene 10m |
| 675 5070P25 | Braided Cord - Polypropylene 25m |
| 675 5070P50 | Braided Cord - Polypropylene 50m |
| 675 5070D04 | Standard heavy duty hand reel |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)



PTFE Mini Dipper

Excellent Chemical Resistance - GMP correct

Ideal for surface sampling small volumes of liquids

- PTFE construction
- Chemically resistant



| Stock No. | Volume | Diameter | Handle |
|-------------|--------|----------|-----------------------------|
| 675 5100A10 | 10mm | 25mm | 155mm long - fixed in place |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)

Disposable Ladle

Efficient surface sampling of liquids



| Stock No. | Volume | Diameter | Pre-Sterilized | Material | Per Box |
|--------------|--------|----------|----------------|----------|---------|
| 675 8053H50 | 50 ml | 48 mm | No | HDPE | 40 |
| 675 8053H50S | 50 ml | 48 mm | Yes | HDPE | 40 |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)

Stainless Ladle

- Satin finish
- 304 or 316 stainless steel
- All stainless steel construction
- Crevice free construction



| Stock No. | Stainless Steel | Bowl Volume | Bowl Diameter | Handle Length |
|-------------|-----------------|-------------|---------------|---------------|
| 675 A326150 | 316 | 150ml | 88mm | 370mm |
| 675 A32450 | 304 | 50ml | 63mm | 275mm |
| 675 A324120 | 304 | 120ml | 82mm | 345mm |
| 675 A324150 | 304 | 150ml | 92mm | 350mm |
| 675 A324200 | 304 | 200ml | 100mm | 365mm |

[QUICK QUOTE LINK](#)

[ONLINE PRICING LINK](#)