


The disposable Powder Sampler is quick and easy to use:

- 1. Push the sampler into the product
- 2. Rotate handle to open up slots
- 3. After leaving the sampler for several seconds to allow powder to enter, rotate handle to close the sampler.
- 4. Remove sampler
- 5. Deposit the sample into a suitable container

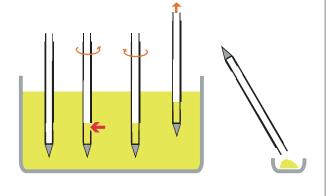

Stock No.	Length	Vol.	Pre-Sterilized	Diameter	Material	Per Box
675 8240H1000	1000mm	100ml	No	40mm	HDPE	20
675 8240H1000S	5 1000mm	100ml	Yes	40mm	HDPE	20

QUOTE REQUEST 💎

BUY ONLINE 💜